



## Assembly Instructions

### 1400 SEAT BACK SET

#30-13-1400

#### Contents of 1400 Seat Back Set

Two Padded Lateral Supports with chest strap

Double Seat Back Support Rods

Two lateral support brackets with cam lock levers

Adjustable Seat Post for AM-12/AM-16 or Angled Seat Post for ProSeries 1412



- 1) Remove installed seat post if this is an upgrade to the seating system. Loosen seat post knob and pull out until seat post can be lifted out.
- 2) Install adjustable seat post provided. Make sure that seat post is pointed toward the front of the trike and it is locked in place. (**Note: Make sure seat post knob is tight before riding**)




- 3) Loosen the black securing knob located in the seat post. This will allow the double seat back support to slide into the bracket.




- 4) Slide the double seat back support rod assembly completely into the black seat support bracket.




- 5) Tighten the black securing knob in seat post.  
(Note: Securely tighten before using)

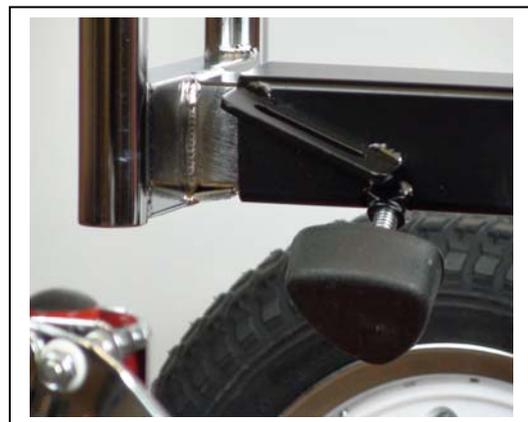

- 6) Make sure the padded seat back can slide freely up and down the double rod supports. Position the padded seat back where desired and tighten the two knobs on the back.

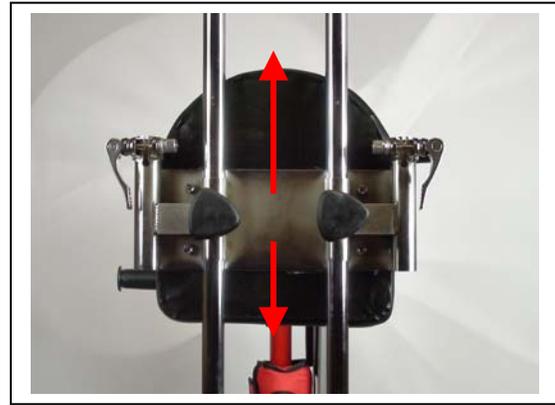

- 7) To position the lateral support brackets, release cam lock lever. The lateral support bracket is now free to move up and down for desired position. Slightly tighten the nut located on the backside of the cam lock lever until you can comfortably close the lever.

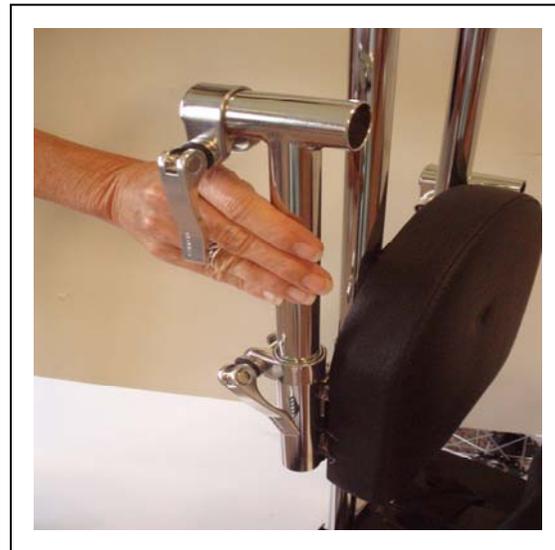

- 8) Lateral supports slide into the lateral support brackets. Adjust lateral supports to desired position and then tighten with cam lock levers.

**NOTE: LATERAL SUPPORTS ARE GROOVED ON BOTTOM AND SLIDE INTO SUPPORTS.**



- 9) Install lap strap from the original seat system. Make sure that hook fastener faces away from rider.



- 10) Install saddle seat onto seat post on seat post assembly. Be sure seat is fully seated onto post. Align with front of tryke and tilt to preferred position. Tighten seat clamp equally on both sides with 9/16 inch wrench.

**Optional Accessory  
Seat Back Spacer Kit  
#30-0502**

Smaller riders may need the seat back position closer to the torso for better support. If so, order the *optional* seat back spacer kit. The kit will allow the seat back to be positioned one and one-half (1 ½) inches closer to the rider.





**Before**



**After**

## Breakdown of components included

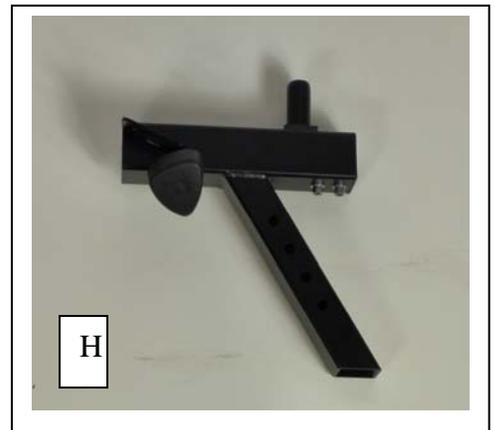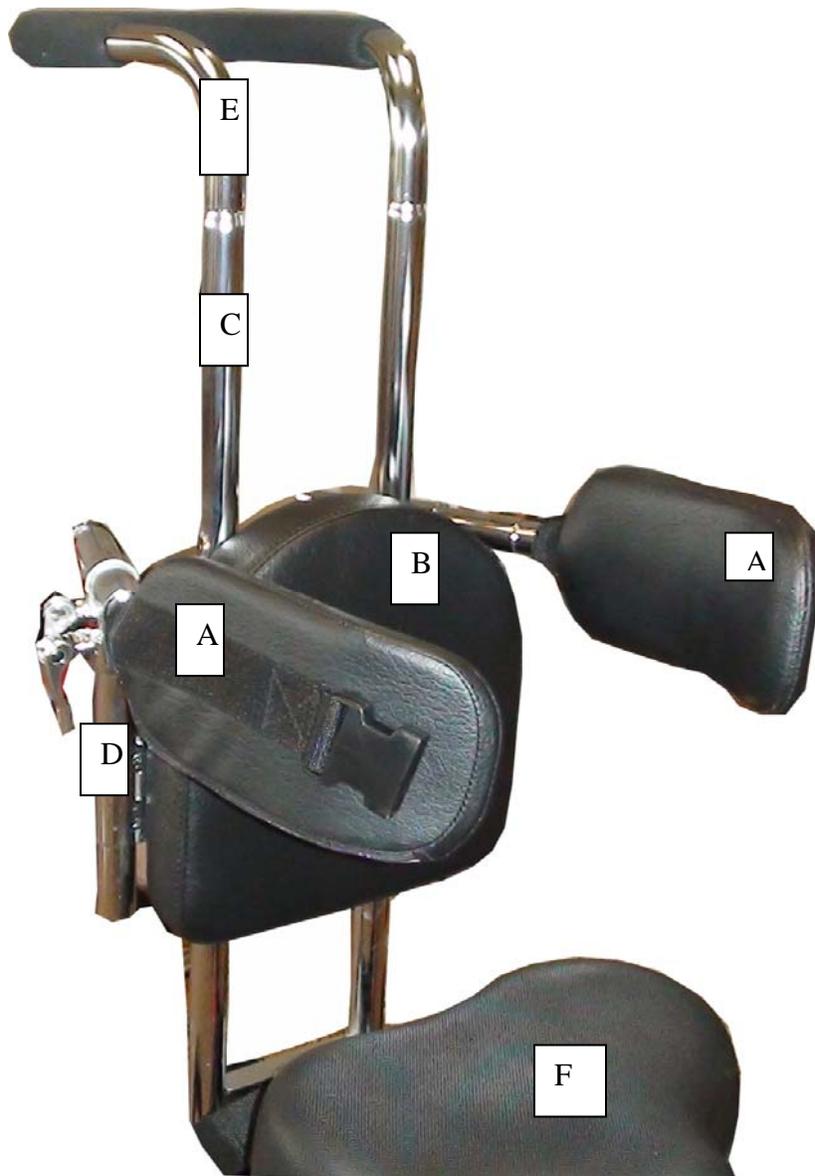

- A) One Pair Padded Lateral Supports (with built-in straps)
- B) Padded Back
- C) Support Rods
- D) Lateral support mounting brackets
- E) Caregiver handle “can be removed to add optional head rest or full padded back”
- F) Saddle Seat. (Saddle seat is not included in #30-0501 1400 seat back set)
- G) Standard Seat Post for AM-12/ AM-16
- H) Angled Seat Post for 1412 (not included unless specified for ProSeries 1412)