

Assembly Instructions

Read these instructions carefully before assembling the product or using it for the first time. They contain important information for your safety and for the correct use and maintenance of this AmTryke® product. Keep the instructions in a safe place for information or for ordering spare parts.

This AmTryke® product has been designed in accordance with the latest standards of safety and subjected to constant quality monitoring. The data gained from this quality process is then incorporated in the development of our products. For this reason, we reserve the right to make changes in design and technology in order to offer our customers products of the highest standard of quality at all times. In spite of this, should you have any ground for complaint, please contact your AmTryke® dealer.

For your safety

- This toy vehicle should be used for its intended purpose only, i.e. one (1) child only. It should not be used in any other way, as this may be dangerous. The manufacturer cannot be held liable for damage or injury caused by incorrect use.
- This product is a toy and is not suitable for use in traffic or on public roads.
- Always keep a watchful eye on your child in order to avoid risks by thoughtless and hasty movements.
- Please care for the fact that the tricycle is only used on safe ground.
- The tricycle is an unlighted play implement. Its use in the dark is dangerous.
- In case of a sudden and/or strong turn of the handlebar the tricycle could be caused to tip. Danger of injuries!
- The tricycle must not be connected with other vehicles, sports or similar devices and not be driven/pulled by them. Never subject your child on the tricycle to speeds exceeding walking pace!
- We recommend that a child using the tricycle should wear a child's bicycle helmet.
- Damaged components may endanger the safety of your child or shorten the lifespan of the product. Inspect it at regular intervals to ensure that none of the components are damaged and that none of the bolts or other joints have worked loose.
- Replace worn or damaged parts immediately and ensure that the product is not used again until this has been done. If any parts require replacement, use only original spares manufactured by KETTLER® and AmTryke®.
- This child's wheeled toy conforms to the stipulated safety regulations Europe standard EN 71-1/3 („Safety regulations for children's toys“). Incorrect repairs or alterations (removal of original parts or addition of other parts etc.) may be a source of danger to the user.

Use of the Toy

- **WARNING:** Adult supervision is required. Never use near steps, slopes, driveways, hills, roadways, alleys or swimming pool areas. Always wear shoes. Never allow more than one rider (without a tandem seat). The use of protective gear and a Safety Helmet is highly recommended.
- It is intended that parents/guardians should seek the qualified opinion of a certified therapist to ensure the AmTryke® therapeutic tricycle is appropriate for use with any child with disabilities.
- Ensure that the toy vehicle is not used until it has been properly assembled and checked by an adult person.
- In case of tricycles without a free-wheel consider the fact that the pedals will rotate even when the tricycle is pushed.
- Do not transport any objects on the tricycle.
- Always adjust the use of the tricycle to the conditions encountered.
- Tricycles do not possess a continuously effective brake system. Therefore, safeguard the tricycle against rolling away when parking it. Never park the tricycle unattended with a child on it!
- Never draw the tricycle behind you as long as a child is seated on it.
- Instruct your child how to use the vehicle correctly and draw attention to possible sources of danger.
- Although dangerous situations are unlikely to occur if the toy is used correctly, children at play should, if possible, be supervised.
- Please bear in mind that there is **always a danger** of children injuring themselves whilst using swings, climbing frames, toy vehicles or similar equipment.

- When using the toy vehicle, the user should always wear strong shoes.
- Instruct your child how to use the vehicle correctly and draw attention to possible sources of danger.
- When children are at play, situations may arise which cannot be predicted by the manufacturer and for which he cannot be held responsible.
- Before use, always check all screws and plug-in connections as well as respective safety devices fit correctly.
- Never use corrosive or abrasive substances to clean the toy. Never use substances which may cause pollution of the environment.
- Lubricate the bearing sleeves of the wheels and the handlebars with a few drops of oil every two or three months.
- For models with push bars: Always ensure that the push bar fits correctly in the adaptor, and locks. Do not hang any weights on the bar! The vehicle may tip over, and your child may be hurt. The push bar is intended for adult use only, do not pass responsibility for children on to other children. For free play time take the push bar off for safety reasons.

Assembly of the toy vehicle

- The toy vehicle must be assembled carefully by an adult person. Study the assembly instructions carefully before beginning work. Assemble the vehicle in the exact sequence shown by the letters in the diagrams. This is essential for the safe function and use of the scooter.
- Please note that there is **always** a danger of injury when working with tools or doing manual work. Therefore please be careful when assembling this machine.
- Ensure that your working area is free of possible sources of danger, for example don't leave any tools lying around. Always dispose packaging material in such a way that it may not cause any danger. There is always a **risk of suffocation** if children play with plastic bags!
- Bolt all the parts together loosely at first, and check that they have been assembled correctly. Tighten the locknuts by hand until resistance is felt, **then use spanner to finally tighten nuts completely against resistance (locking device)**. Then check that all screw connections have been tightened firmly. **Attention:** once locknuts have been unscrewed they no longer function correctly (the **locking device is destroyed**), and must be replaced.
- The grips can be pushed onto the handlebars more easily if they are first soaked in warm water.
- Please keep original packaging of this article, so that it may be used for transport at a later date, if necessary.
- For technical reasons, we reserve the right to carry out preliminary assembly work (e.g. addition of tubing plugs).
- When ordering spare parts, always state the full article number, spare-partnumber, the quantity required and the inspection number stamped on the back.
Example order: Art. no. 08844-999 / spare-part no. 12802003 / 1 piece / inspection no. ...
Important: spare part prices do not include fastening material; if fastening material (bolts, nuts, washers etc.) is required, this should be clearly stated on the order by adding the words „with fastening material“.

Warranty

3 Year Residential Warranty and a 1 Year Institutional Warranty.

The warranty begins with the date of purchase by the consumer. Please address all warranty questions to AmTryke® LLC, 800-838-1845.

For parts or service ONLY call 800-838-1845 or send e-mail to:

World Wide Web address: www.ambucs.org

Email address: ambucs@ambucs



The manufacturer reserves the right to make changes reflecting product or other technical advances. The content of these instructions is accurate at the date of printing.

Checklist


	1	✓
	1	
	1/1	
	1	
	1	
	2	
	1	
	1	
	1	
	1	

		1	
		1	
		1	
	M12	1	
	Ø 6 mm	4	
	M6x65mm	2	
	M5x30mm	2	
	4,2x16 mm	5	
	4x16 mm	2	
	4,8x16 mm	1	
		1	
		1	
		1	

1

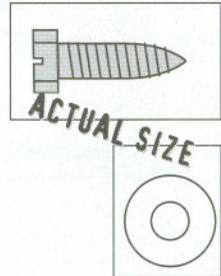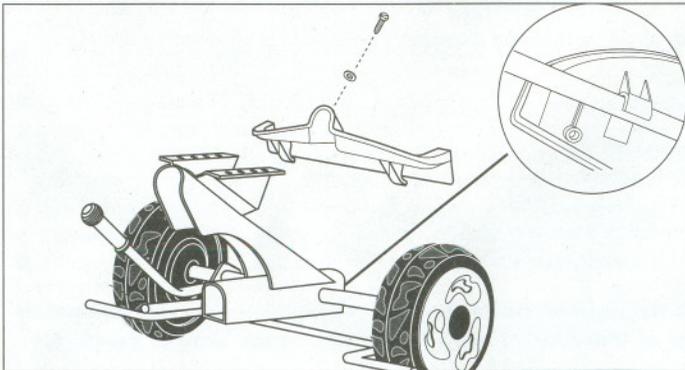


2

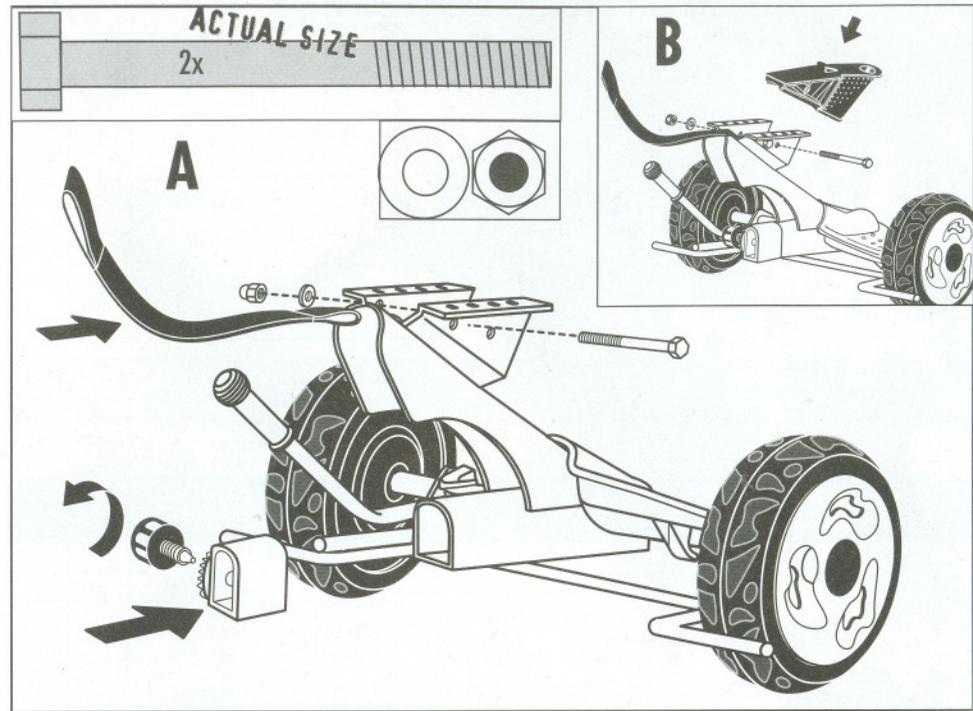


ALTERNATE:



Without safety belt and / or without adaptor the screw holes have to be closed with the plugs.

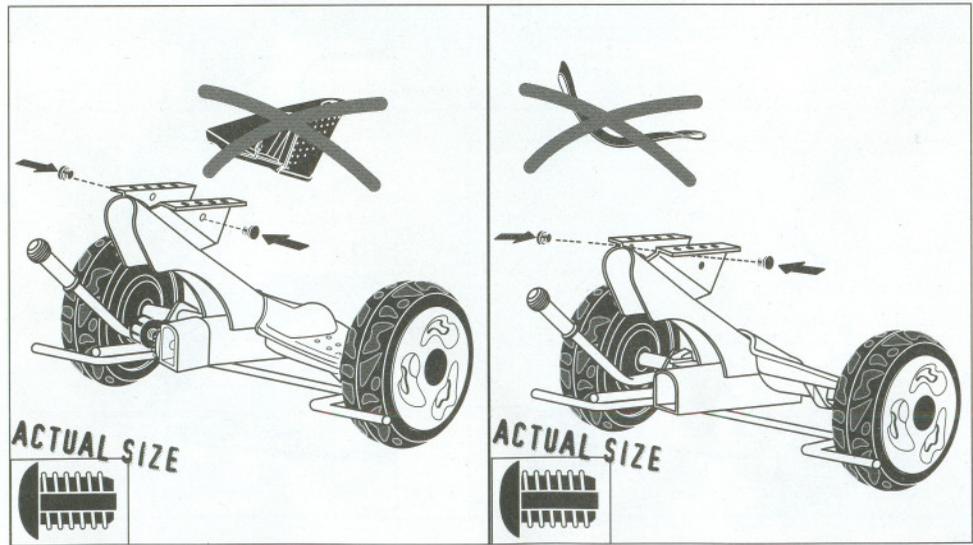

3


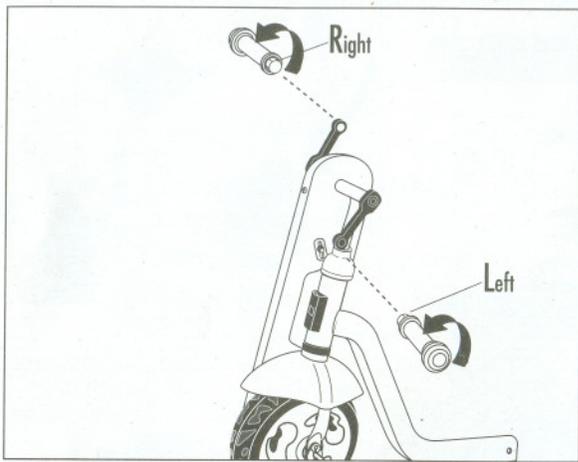


4

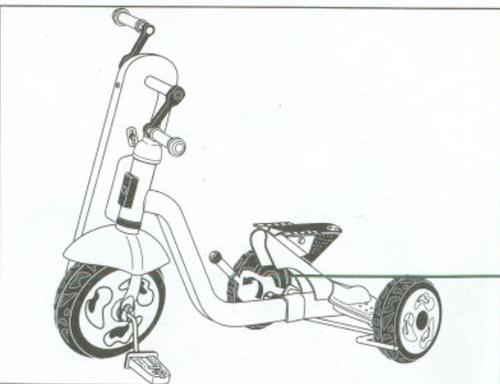
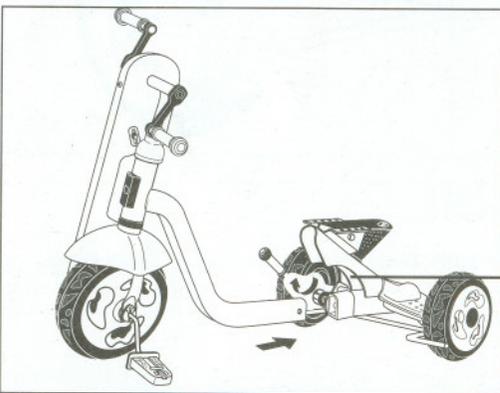


6





5

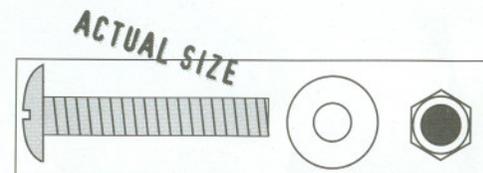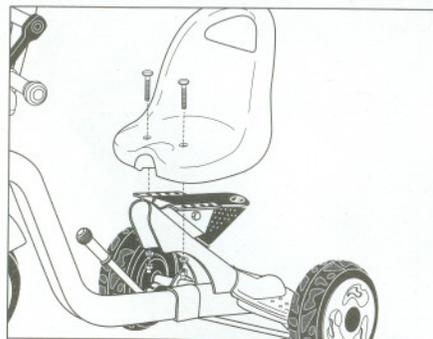

6

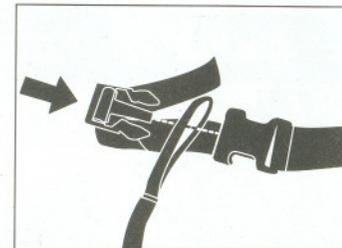
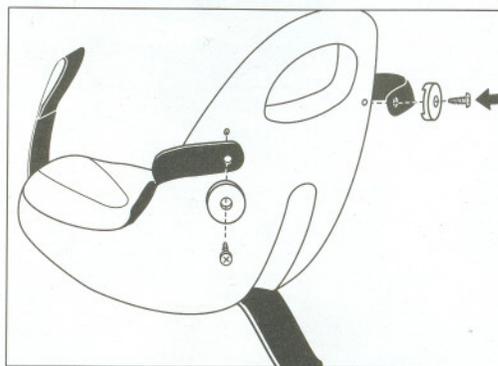


7

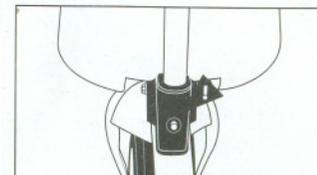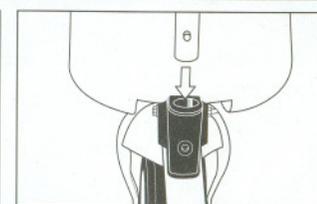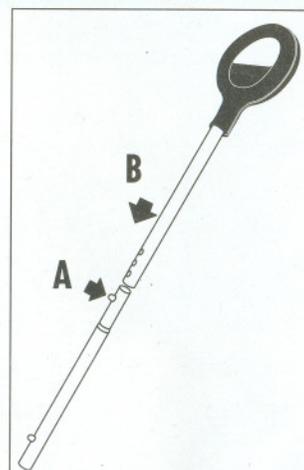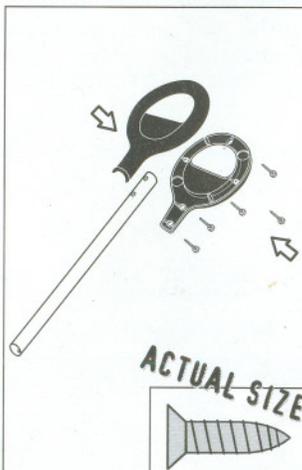


7


Handlebar locking/unlocking

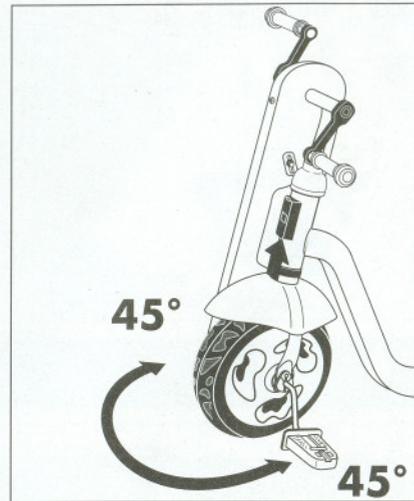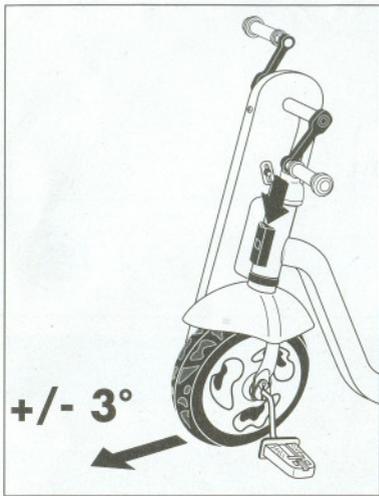

Tightening/loosening the chain

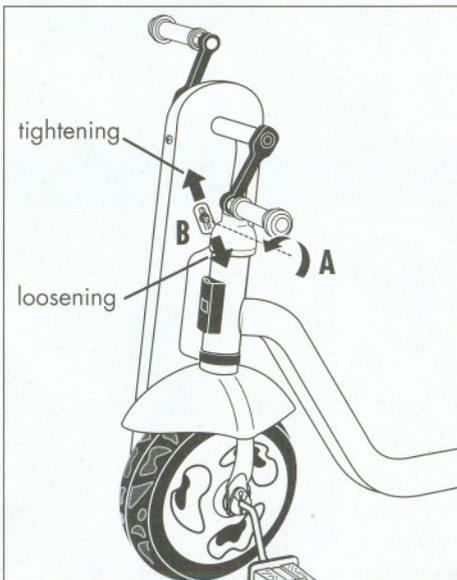

Spare parts drawing

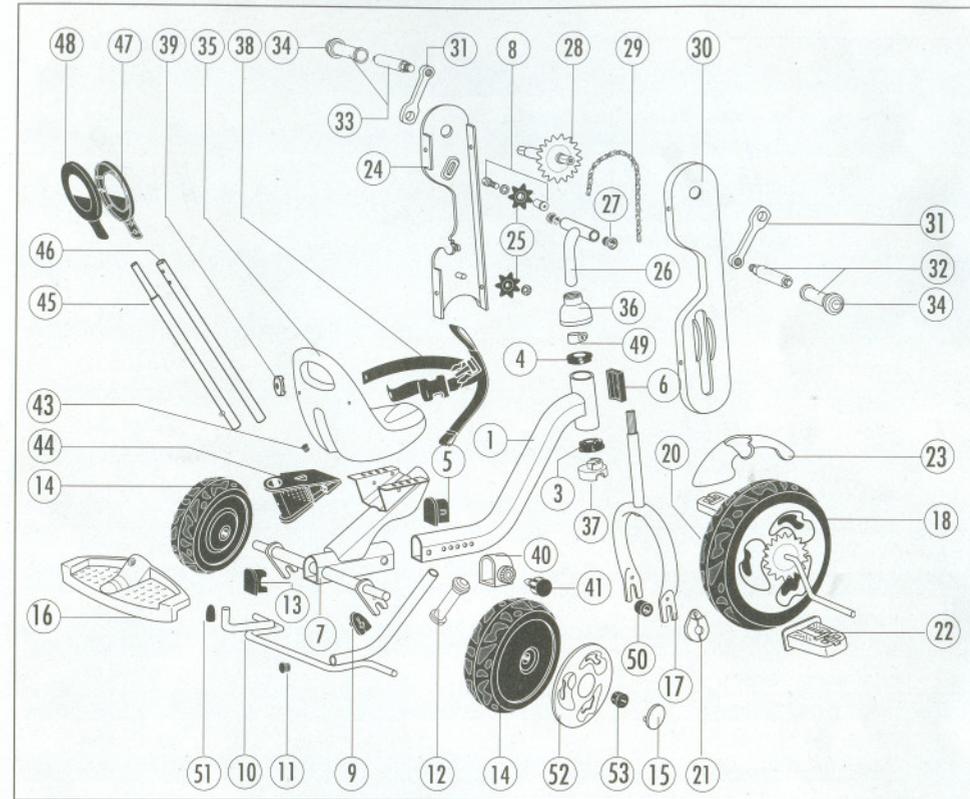

When ordering spare parts, always state the full article number, spare part number, the quantity required and the inspection number stamped on the back.

Example order: Art. no. 08844-999/spare part no. 10100169 /1 piece/inspection no.

Important: spare part prices do not include fastening material; if fastening material (bolts, nuts, washers etc.) is required, this should be clearly stated on the order by adding the words „with fastening material“.

3 Year Residential Warranty and a 1 Year Institutional Warranty.

The warranty begins with the date of purchase by the consumer. Please address all warranty questions to AmTryke® LLC, 800-838-1845.

For parts or service ONLY call 800-838-1845 or send e-mail to:

World Wide Web address: www.ambucs.org

Email address: ambucs@ambucs

pos.	description	qty.	part no.
1	front frame complete	1	91611040
2	locking pin for parental control (cpl., w/o fig. 2968)	1	91640021
3	bottom headtube bushing for ø42mm (2969)	1	70121560
4	headtube bushing (2046)	1	70122051
5	front frame plug 35x41 mm (2995)	1	70128425
6	parental control locking mechanism (2970)	1	70121770
7	rear frame	1	91611006
8	nut + bolt for chain tension	1	10321002
9	axle bearing	2	10116029
10	brake lever	1	91640032
11	plug	1	10100006
12	grip for brake lever (2517)	1	70121571
13	frame plug (3535)	1	70121432
14	rear wheel	2	15170087
15	plug for wheel sceen (2042)	2	70121823
16	foot plate	1	70122020
17	fork	1	91620293
18	front wheel	1	15170090
20	hub cap left (3918) w/o fig.	1	70121452
21	front wheel bearing (3920)	2	70121453
22	pedal for Ø 8 mm (1848)	2	70121541
23	fender (3921)	1	70121454
24	chainguard, lower part	1	91670109
25	sprocket wheel	2	33101101
26	handlebar	1	91650267
27	bushing (3919)	2	70121455
28	axle 13x150	1	97201894
29	chain	1	42208084
30	chainguard, upper part (3914)	1	70121456
31	pedal arm (left + right)	1	33001012
32	handlebar axle with grip, right	1	91630083
33	handlebar axle with grip, left	1	91630084
34	grip	2	10111023
35	molded seat (3988)	1	70120244
36	covering (3539)	1	70121357
37	handlebar limiter (3536)	1	70121397
38	safty belt	1	91680284
39	fastening for safty belt (3157)	2	70121640
40	adjustement sleeve (3064)	1	70121435
41	quick release knob M12	1	91670058
42	hardware bag (w/o fig.)	1	91680358
43	plug 8,5x6,5x12 mm	4	10100169

pos.	description	qty.	part no.
44	adaptor for push bar (3493)	1	70125638
45	pushbar, lower part	1	97100400
46	pushbar, upper part	1	97100399
47	grip for pushbar, upper part (3989)	1	70121575
48	grip for pushbar, lower part (3379)	1	70121574
49	c clamp	1	97200347
50	PE distance roll	1	10108160
51	lacquerd cap 7,5x25 mm	1	10101105
52	hub cab for running wheel	2	70121845
53	wheel bearing	4	70121470

AmTryke®, LLC is a wholly owned subsidiary of National AMBUCS™, Inc. AMBUCS™ is dedicated to creating independence for people with disabilities, which includes providing AmTryke® to children with special needs.

Check AMBUCS™ out on-line at www.ambucs.org!